



# Etac Clean 24"

For the independent user – easy to use and propel.

## Easy manoeuvring

The rear wheels on Clean 24" are positioned to maximise the manoeuvrability which makes the chair easy to propel without compromising the security of the user.



## Easy transfers

Side transfers can be performed safely as there is sufficient seat space in front of the wheels. Furthermore the brake handles are positioned just underneath the seat surface so they do not obstruct transfer but are still within easy reach.



*There is enough seat space in front of the wheels for side transfer.*

*An Etac Clean with 55 cm (21¾") seat height can be remodelled to a Clean 24" by assembling the accessory "Rear Wheel Adapter".*

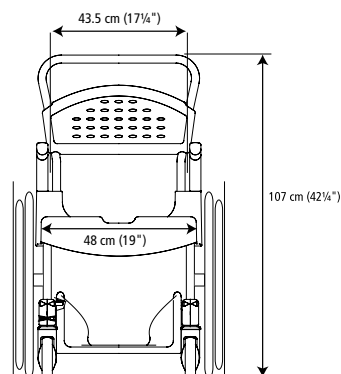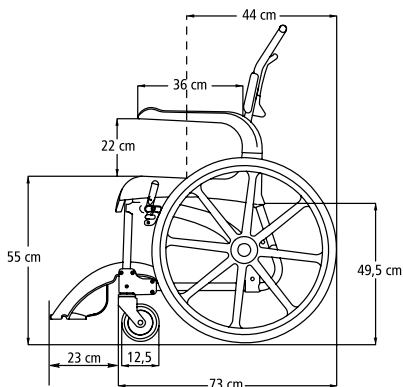

## Accessories

See our wide range of accessories "Accessories for all Clean models" on page 24

### Description

Clean 55 cm (21¾) with 24" rear wheels, green 80209216  
 Clean 55 cm (21¾) with 24" rear wheels, grey 80209276

### Material

Rear wheel adapter: aluminum  
 Rear Wheels: polyamide, solid tires  
 Other material: see page 18

### Weight

20 kg (44.1 lbs)



Clean the product with a solvent-free cleaning agent (pH level 6–9), or with a 70% disinfectant solution. Can be decontaminated at max. 85°C / 185°F if so required.